



# 9G!7 @ \$

7 C & Gi f [ ] W U ~ U g Y f ' g m g h Y a



OPERATING THEATRE



# EG!7 @ \$

## CO<sub>2</sub> Surgical laser system



**Cod. ES-CL30**


- Sealed-off laser technology
- Cw. pulse and ultra pulse operation modes
- Touch switch, microprocessor message display
- 7-joint spring-balance articulated arm
- 635nm diode laser pilot beam
- Safety protecting or with automatic alarm
- Control panel:



ES-CL30 CO<sub>2</sub> Surgical laser System is equipped with a reasonable construction, stable performances, safety and reliability. The CO<sub>2</sub> surgical laser therapy apparatus is suitable for different departments as dermatological, ENT, gynaecological and general surgical for cut, cautery and vaporization treatment. This machine is used by different levels of hospitals and medical centers.

ES-CL30 CO<sub>2</sub> Surgical Laser has got an ultra pulse function: the laser can reduce hot damage around the target organ and the output time is extremely short (1 ms). This device has a blowing function which around the target organ can be cooled at pulse interval time. Adjustable interval time of ultra pulse that controls the depth of laser act.

ES-CL30 CO<sub>2</sub> Laser Surgical System is also suitable for pixel beauty treatment, such as removal pigment, wrinkle, blemish and rejuvenation.

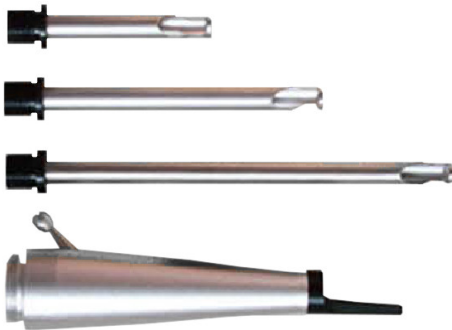

### Standard enclosure head pieces

F=100 mm head x1  
Lengthening head x3  
Pixel head x1



### Optional head for operation



ESSE3 srl, Via Garibaldi 30  
14022 Castelnuovo D.B. (AT)  
Tel +39 011 99 27 706  
Fax +39 011 99 27 506  
e-mail [esse3@chierinet.it](mailto:esse3@chierinet.it)  
web: [www.esse3-medical.com](http://www.esse3-medical.com)

