

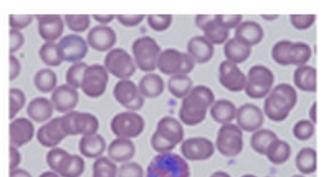
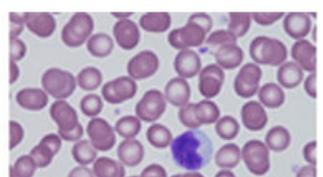
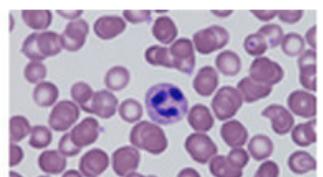
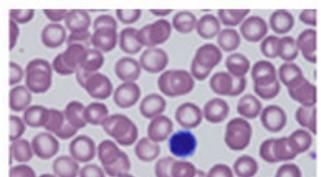
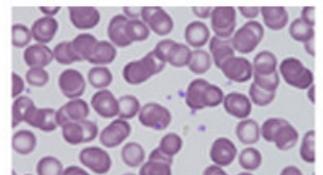
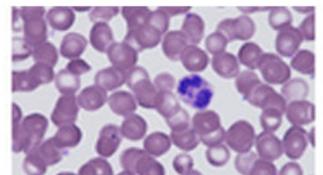


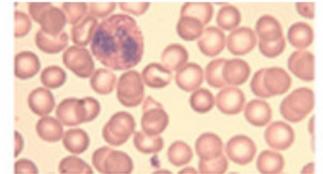

AUTOMATED SLIDE STAINER

Cod. QuickSlide Plus II

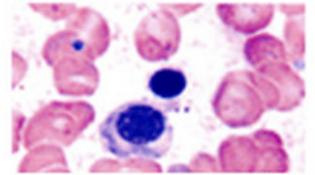

QuickSlide Plus II AUTOMATED SLIDE STAINER

The QuickSlide Plus II is one of the smallest automated slide stainers available and it can stain peripheral blood specimens in a time of 45 sec. The stainer features a microprocessor to help increasing speed and a program that enables it to stain bone marrow samples. The stainer helps to eliminate cross contamination associated with the dip method and offers consistent staining from user to user. Other features include an LCD screen operator, soft-touch color-coded keys, adjustable stain/buffer timing and fill level, cycles remaining display and reagent monitoring capabilities. Housed in an injection-molded case with a built-in reagent platform, the stainer also includes a high-frequency transducer speaker with audible sound that requires less technician time. The unit's stainless steel reagent cannulas are oxidation-resistant to help increasing system durability.

OPTIONAL ACCESSORIES

QS-29 Reagent Pack

Reagent kit sufficient for 125 tests

QS-2 Retubing Kit

The retubing kit includes tubing, elbow fittings, pump windings, stainless steel cannulas, one-way valves and wire ties for a complete front and back replacement.

TECHNICAL SPECIFICATIONS

Cod. QuickSlide Plus II

Automated Slide Stainer

Displays	"Cycles remaining" displayed; computer monitors reagent levels
Required accessories	Requires no pre-fixing of blood smears; bone marrow smears require pre-fixing with methanol
Adjustability mechanism	Independantly adjustable stain & buffer times to customize intensity
Multicolor	Color coded tubing for easy reagent replacement
Temperature, ambient operation	16°-32° C (60°-90° F)
Humidity	0 to 95% without condensation
Operator adjustments	Stain/Buffer Times Adjustable Fill Levels On & Off Audible Sound Control
Operator entry	Moisture Proof, Soft Touch Keypad Stain Resistant Injected Molded Cabinet
Width	30.48 cm
Height	27.9 cm
Depth	27.94 cm
Weight	4.19 Kg
Power requirements	Input 100-240 50/60HZ

Cod. QS-29









Reagent Pack

Estimated Number of Slides Stainable	125 slides
Use With	QuickSlide Plus II
pH Range	pH of the buffer is right around neutral of 7.0 (between 6.8 and 7.2)
Average Shelf Life	Unopened; 1 year shelf life prior to opening
Average Shelf Life	When opened is good for 3-7 months; can be used as long as the QC is acceptable

Cod. QS-2

Retubing Kit

Use With	QuickSlide Plus II
Suggested Use	Recommend tubing replacement every 6-8 months of use

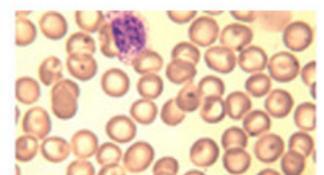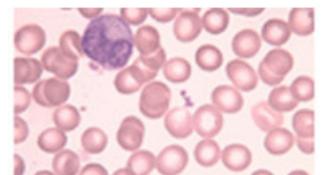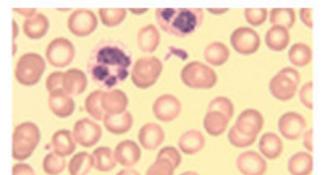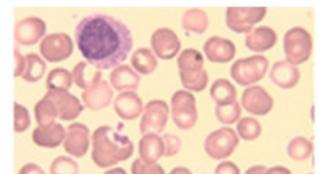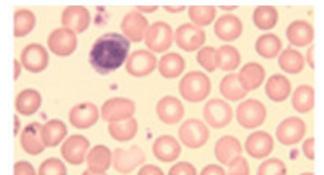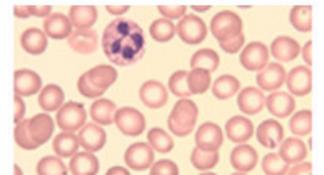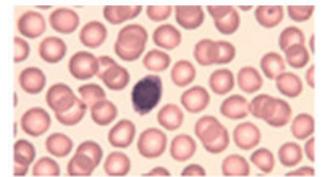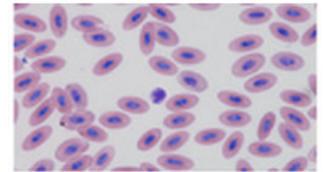

ESSE3 srl, Via Garibaldi 30
14022 Castelnuovo D.B. (AT)
Tel +39 011 99 27 706
Fax +39 011 99 27 506
e-mail esse3@chierinet.it
web: www.esse3-medical.com

