



ES-NT1D



Handheld Capnograph / Pulse Oximetry Monitor



Cod. **ES-NT1D**

- Durable, Compact and Lightweight
- Capnograph with Trends
- SpO2 with Waveforms
- Data Storage for up to 100 patients
72-Hour for Each Patient
- Wireless Data Transmission
- Suitable for Adult, Pediatric
and Neonate





Portable

Compact & Ergonomic Design




21st Century CO2 Technology

CO2 sensor and
LoFlo Sidestream System



All-in-One Display

Data + Waveform + Trend

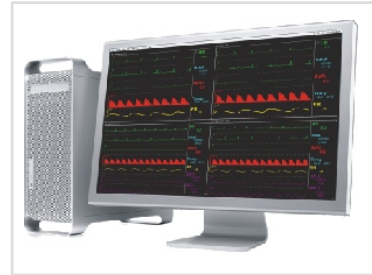

Data Management

PC Software with Wireless Data Transmission
Data Analysis and Report Printing



A Variety of Applications

Emergency Rescue, Intensive Care
During Surgery, Resuscitation and
Patient Transportation



For All Patient Types

Adult, Pediatric and Neonate



Usage Environment

- Emergency Medical Services (EMS) in the field or during transport
- Outpatient or Ambulatory Surgery centres; special procedures area (e.g. cardiac catheterization lab, endoscopy)
- General medical/surgical hospital ward
- ICU, Emergency Department
- Hospital-based or free-standing sleep laboratory

Clinical Applications

- Airway management for all intubated patients
- Procedural or conscious sedation-adequacy of ventilation
- Patient safety during patient-controlled analgesia (PCA) or continuous narcotic administration
- Cardiopulmonary resuscitation-- confirm endotracheal tube placement, determine effectiveness of chest compressions (CPR) and detect Return of Spontaneous Circulation (ROSC)
- Sleep Studies

Why Capnography is A Valuable Tool

- For EMS Transport:
Capnography is a valuable tool during emergency transport of both intubated and non-intubated patients for proper assessment of the patient's ventilatory status.
- For Conscious Sedation:
When performing procedural sedation, ensuring patient safety and adequate ventilation is essential
- For Cardiopulmonary Resuscitation:
Capnography is a valuable tool during cardiopulmonary resuscitation (CPR) of intubated patients
- For Pain Management:
The use of capnography is becoming more widespread for patients receiving opiates for acute pain management
- For Sleep Laboratories:
When conducting sleep studies, it is important to accurately and consistently measure exhaled CO₂ levels in order to reliably assess the quality of ventilation during sleep

Available SpO₂ and CO₂ sensors



Adult / Pediatric Finger



Adult / Pediatric Soft-Finger



Single Patient Disposable



Mainstream CO₂



Sidestream CO₂



Disposable Cannulas

Technical Specifications

SpO₂:

Measurement Range: 0 ~ 100%
Accuracy: +2% during 70%~100%
0%~69% unspecified

Pulse Rate:

Measurement Range: 30 bpm ~ 250 bpm
Accuracy: 1 bpm or $\pm 2\%$
whichever is greater

EtCO₂:

Measurement Range: 0~150mmHg
Resolution: 0.1mmHg (0~69)mmHg
0.25mmHg (70~150)mmHg
Accuracy: ± 2 MmHg (0~40)mmHg
 $\pm 5\%$ (41~70)mmHg
 $\pm 8\%$ (71~100)mmHg
 $\pm 10\%$ (101~150)mmHg

Respiration Rate:

Measurement Range: 0~150bpm
Accuracy: ± 1 bpm

Alarm:

Three levels of visual, audio alarms

Data Transmission:

2.4GHz wireless
USB to PC
 ≤ 10 m without obstruction

Power Requirements:

DC: 9V 700mA
4x2400mA rechargeable batteries
Battery Capacity: ≥ 12 hours (SpO₂ only)
Battery Capacity: ≥ 5 hours (SpO₂ + CO₂)

Environment:


Operating Temperature: 0°C~50°C
Humidity: $\leq 95\%$
Altitude: -390m~5,000m
Transport/Storage Temperature: -20°C~70°C
Humidity: $\leq 95\%$

Physical Characteristics:

Dimensions: 73mm (W) x 127mm (H) x 23mm (D)
Maximum Weight: 500g



Display Options

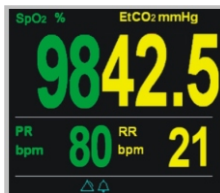

Large Font/Digits



Dural Waveforms



Historical Trend



Trend Chart

Configurations

- ES-NT1D-B Handheld Mainstream CO₂ Monitor
- ES-NT1D-C Handheld Sidestream CO₂ Monitor
- ES-NT1D-D Handheld SpO₂ & Mainstream CO₂ Monitor
- ES-NT1D-E Handheld SpO₂ & Sidestream CO₂ Monitor



ESSE3 srl, Via Garibaldi 30
14022 Castelnuovo D.B. (AT)
Tel +39 011 99 27 706
Fax +39 011 99 27 506
e-mail esse3@chierinet.it
web: www.esse3-medical.com

