



EST-90

Universal Hydraulic
Operating Table



OPERATING THEATRE





EST-90

Universal Hydraulic Operating Table



Designed to meet different surgical requirements, the EST-90 is a multi-function operation table with HEAD-END control system. Its design is suitable to every requirement.

- Stable base construction with four castors for easy movement.
- Base brake lock/unlock by footswitch.
- Lateral tilt and Trendelenburg adjustable by hand cranks.
- Base, side rails, accessories & clamps made of #304 stainless steel
- Table top elevation controlled by hydraulic system
- Table top was divided into 5 sections including head plate, back plate, seat plate and leg plates
- 360° table top rotation
- Kidney Surgery position
Table top can be adjusted to the Flex/Reflex positions for kidney surgery

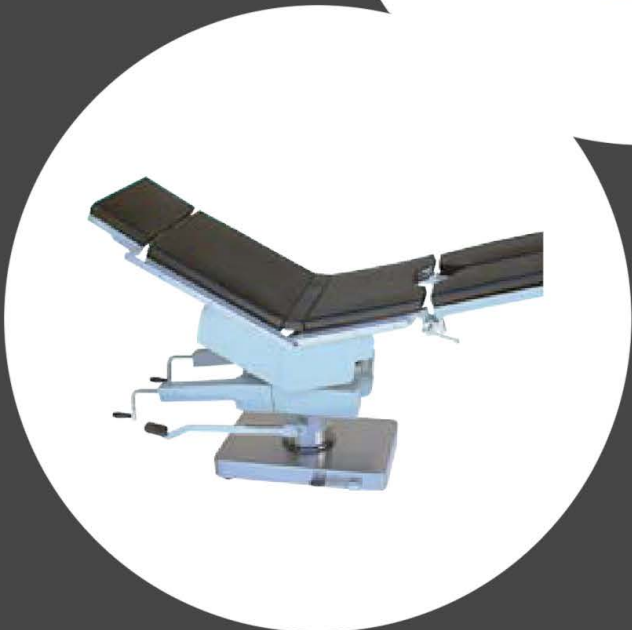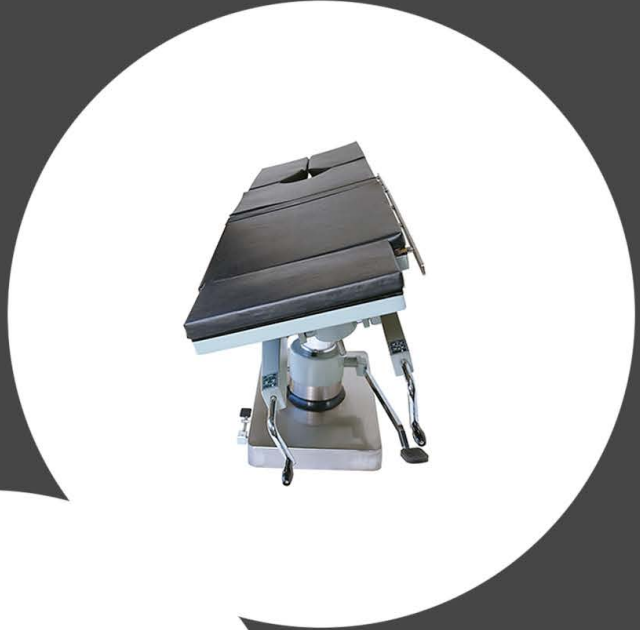



Technical Features

Table-top length	200cm
Table-Top width	50cm
Table-top elevation	68 - 100cm
Left/right lateral tilt	25°/25°
Trendelenburg/Anti-Trendelenburg	30°/22°
Back section up/down	+80° / -20°
Head section up/down	+30° / -90°
Leg section up/down	+20° / -90°
Seat section up/down	+20° / -20°

Standard Accessories

Screen frame	1 set
Shoulder supporter	1 pair
Body support	1 pair
Armrests	1 pair
Knee crutches	1 pair
Mattress	1 set

Optional Accessories*

Restrain strap	1 pc
Hand grips	1 pair
U-shaped headrest	1 set
Raised armrest	1 set
X-ray radiograph table-top	1 set
Contamination bucket	1 set

*For the full optional accessories list, please check the dedicated catalog



ESSE3 srl, Via Garibaldi 30
14022 Castelnuovo D.B. (AT)
Tel +39 011 99 27 706
Fax +39 011 99 27 506
e-mail: esse3@chierinet.it
web: www.esse3-medical.com

